

MORTON


MN 120 | **Coaxial Circuit**
Coaxial Devre

Features:

1. Compact design for ventilator and anesthesia equipment.
2. The circuit can be transferred to the ventilator in POR or ICU during the anesthesia process.
3. Coaxial tubing utilizes the temperature of exhaled air to increase and maintain the temperature and moisture level of the air to be inhaled in the interior tubing (blue).
4. Materials conform to ISO-10993 bio-compatibility requirements.
5. Single-patient use to minimize the risk of cross-contamination.
6. Optional bacterial filter is available upon request.

Specification:

EVA corrugate tubing or 150cm with patient connector 22/15mm.
T Piece EVA foldable tubing 22F

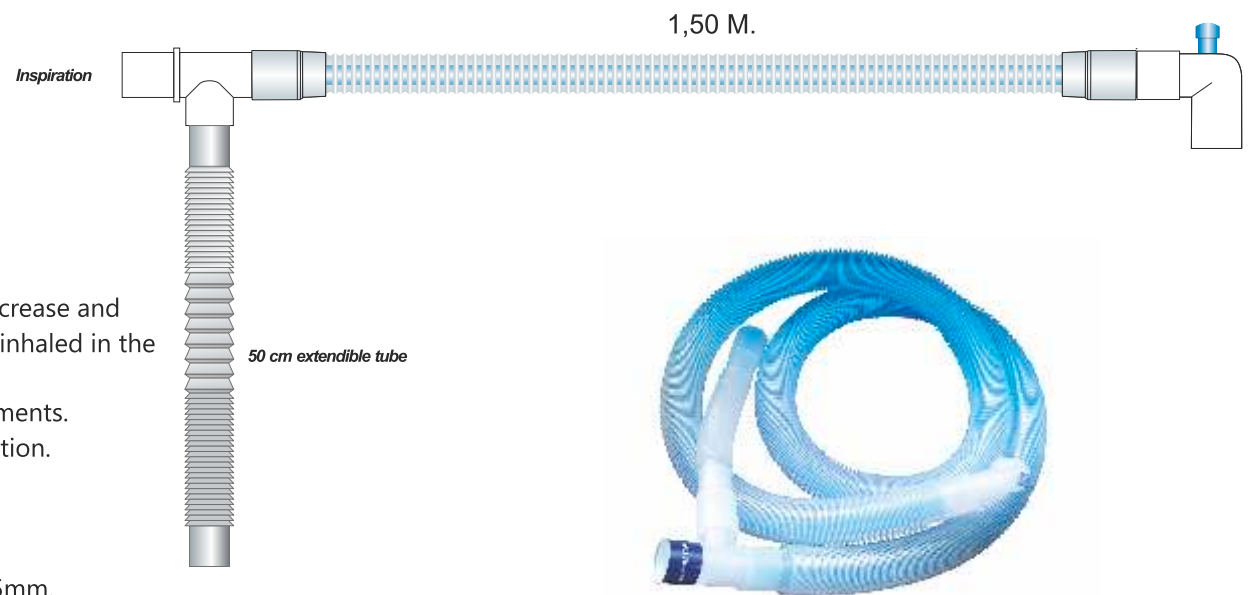

Also MN 208 Gas Sample Line & MN 124 Disposable Rebreathingbags are available upon request.